




College Park designs and manufactures a full line of anatomically correct, custom built prosthetic solutions. College Park is ISO 13485 certified with an extensive focus on quality testing, and proudly sources and produces all products entirely in the USA. With precise engineering and quality manufacturing, College Park is committed to innovating human locomotion solutions for users all over the world.



 **AURA**™
polycentric knee




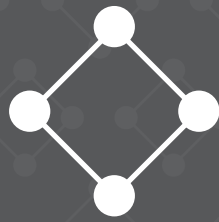

At College Park, our ideal Aura is comforting, intuitive, and well-grounded. The new Aura polycentric knee was designed to offer stability and relief to low-impact users. Featuring a unique locking mechanism and simple 4mm adjustments, the Aura goes beyond mere physical support. It can open doors for mobility goals and newfound peace of mind.

UNIQUE MANUAL LOCK

Users can easily lock and unlock through the simple push of a button. This unique method may improve the user experience and enable usage in water.

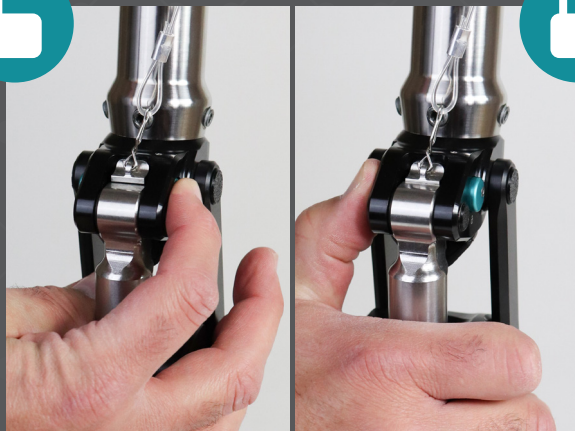

PUSH BUTTON LOCK

REMOTE LOCK



GEOMETRIC STABILITY

The Aura was created to offer inherent stability. With its poly-axial design, the Aura's Instant Center of Rotation (ICR) stays posterior to the load line in stance to provide users with enhanced stability—while still allowing smooth and effortless flexion into swing. When the Aura knee flexes, it shortens, which may increase toe clearance and safety when walking.

BETTER TOGETHER

Pair the Aura with any K2 foot from College Park to maximize product performance! Whether you're seeking tri-axial comfort, hydraulic movement, or enhanced stability, we have solutions for you.




REFINED ADJUSTMENTS

Like all College Park knees, adjustments are quick and easy. Clinicians can modify stance flexion, extension assist, and more with a 4mm hex key.



WATERPROOF

Designed and tested for corrosion resistance, the Aura has been approved for submersion in the ocean and chlorinated pools. It was built to perform well and maintain high quality in different wet environments.

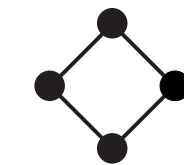

TECHNICAL SPECIFICATIONS

TYPE <i>Polycentric</i>	WEIGHT 665 g
FLEXION RANGE 0-135°	WARRANTY 2 years
OVERALL HEIGHT 6.32" (161 mm)	L CODES* L5611, L5845, L5850, L5925
WEIGHT LIMIT 300 lbs (136 kg)	
DOMES TO KNEE CENTER .65" (17 mm)	IMPACT LEVEL
WATER CAPABILITY Waterproof	

*The listing of codes with these products should not be construed as a guarantee for coverage or payment. Ultimate responsibility for the coding of services/products rests with the individual practitioner.

